

# Dental Unit Line

CARE 11  
Top mounted



**Clamshell Type Film Viewer**  
When film viewer is closed, the space of reaching medicine is more open

**Automatic Check on the Rack**  
Micro switch with automatic fault checking function  
Durable silicon tubes for all instruments

## Clean and Sanitary

Silicon pad on the instrument tray—  
Protection from chemical corrosion  
Autoclavable to avoid cross-infection



**Runyes®**  
FOCUS ON DENTAL

Ningbo Runyes Medical Instrument Co.,Ltd.

Add:No. 456, Tonghui Road, Jiangbei Investment & Pioneering Park C, Ningbo  
Tel:+86-574-27709922 Fax:+86-574-27709923  
Http://www.runyes.com E-mail:runyes@runyes.com

With automatic cleaning,  
Spend your day leisurely;  
With water gurgling,  
Enjoy your mind traveling!

# Dental Unit Line

CARE 11



Perfect your clinic!

**Runyes®**  
FOCUS ON DENTAL



Adjustable Lamp  
Light and flexible  
Model:D2

Nurse Station  
Water-air exchange function makes  
operation more convenient.

Comfortable and Convenient  
90° rotational basin  
---comfortable and convenient

Linkage Structure  
Chair position can be adjusted;  
Chair-back will move as the cushion moves.

Safety  
With protection device to make sure the  
operation safety

Stainless Steel Foot Control  
Stepless speed adjustment  
Available control of mouth wash and basin water

Foot Pedal  
Durable PU material  
Available control of chair up-down movement

Adjustable Mouth Wash Temperature  
---more care for the patients

## Dental Unit Line CARE 11 Down mounted



FARO LED lamp  
---available for option  
---Model:D4




LED lamp  
---available for option  
---Model:D3



Creative Rack  
Detachable parts of the rack  
Adjustable space from the board.

Silicon Tube  
Selected material  
Durable for years

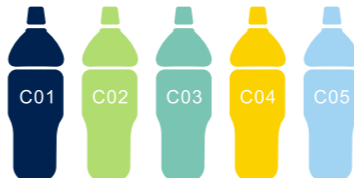

Detachable filter for easy cleaning  
Dirt separator---avoid jam



Top location of storage bottle mouth, easy  
for refilling. Meet the demands of water-  
break and pipeline sterilization



Double joint headrest: it can bring patients  
a comfortable position and can be used for  
special patients.



Luxury Doctor's Chair ( For standard  
configuration)



Luxury Nurse's Chair ( For Option )

## Dental Unit Standard Configuration and Features

Dental unit standard configuration:		Features:
Handpiece tube	3sets	Imported pipe
3-way syringe	2sets	With constant temperature water supply system
Water heating system	1set	Temperature range: 30℃ ~ 45℃ (Optional)
Suction system	1set	Time delay function
Strong suction	1set	Time delay function
Oral halogen lamp	1set	8000-15000Lux
Film viewer	1set	Special clamshell type film viewer
Rotatable side-box	1set	Equipped with a revolving ceramic basin
Cleaning system	1set	Top location of water mouth, easy for refilling
Air lock type arm	1set	Electronic controlled air lock switch
Footswitch	1set	New type multi-functional footswitch
Dentist chair	1piece	Human body engineering design, beautiful and practical
Water and electricity general switch	1 set	Easy for control
Optional:	Italy 3-way syringe/ high and low speed hand piece/ oral camera/ ultrasonic scaler/ LED curing light / fiber optic/LED lamp/monitor holder	

## Dental Unit Standard Configuration and Characteristics

Dental Unit Standard Configuration:	Characteristics:
Power-driven dental chair	1. Comfortable dental chair of linkage system and human body engineering design 2. Safety protection design
Cushion and backrest	Seamless PU material of special craft---clean, soft and satisfying
Digital computer program setting	Three chair position user-defined memory
Integrated PU base	Perfect combination of luxury and steady.

Technology Parameter:	Packing/Weight	Packing Dimension(Carton)
Power supply AC 220V 50Hz/60Hz	Dental chair	1450x900x935mm
Water pressure 0.2 MPa~0.4MPa	Dental unit	1610x1130x1340mm
Air pressure 0.55MPa~0.6MPa	Dental unit	Net weight: 226 KG Gross weight: 336 KG

## Clinic layout

